

**HMG
DCS**

DEFIBRILLATOR MONITOR

Provide a full range of functions to meet various life support needs

With Defibrillation, Pacing, Monitoring and AED mode, SA is not only suitable for pre-hospital first aid, but also applicable in-hospital usage.

Defibrillation

Manual defibrillation modes include synchronous cardioversion and asynchronous defibrillation.

Pacing

Having on-demand pacing and fixed pacing mode, for patients with cardiac arrest and acute severe slow arrhythmia, in vitro non-invasive pacing mode is rapid, easy to master, time-saving and improve recovery success rate.

Monitor

A variety of monitoring functions including ETCOR, IBP, and 12-lead ECG are available for continuous monitoring of patient vital signs.

AED

The model applies patented analysis algorithm and automated analysis as well as convenient setting to guide clinical emergency personnel in providing defibrillation and basic life support.



**HMG
DCS**

The Most Convenient and Efficient Life Rescue

As the most important part of CPR, time is the key for defibrillator monitor. Therefore, abandons complex operation and improve convenience and efficiency for rescue.

1 for 2

electrode can be divided into large and small electrode paddles which can defibrillate for both adult and child respectively.

1 button

Selecting 12-lead monitoring interface.

3 Steps

fulfill defibrillating operation
(Energy Selection, Charging, Discharging).

Higher and Greater for Better Rescue

Up to **360J** energy selection, for patients with high defibrillation thresholds such as myocardial infarction, obesity, high impedance, higher energy selection indicates greater defibrillation success rate. Greater Impedance Range from **20– 250Ω**, suitable for a great range of patients. Higher efficiency with Advanced Biphasic Truncated Exponential (BTE) waveform technology and automatic impedance compensation.

Standby as you need, review as you request

Reliable quality has always been the pursuit of COMEN R&D personnel. As a first-aid device, defibrillator is often utilized in extreme environments, whereby its reliability is highlighted in various circumstances.

Information-based Network Solution

Provides a comprehensive networking solution, allowing users to access patient information and improve work efficiency when needed.



Lower energy




Lighter damage



Better defibrillation



Anti-shock and anti-fall, sturdy and durable. With IP44 ingress protection, S8 avoids the penetration of liquid and can be applied in complex outdoor environment.

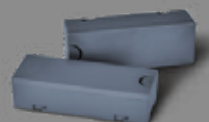

Large-capacity lithium battery, support fast charging and more than 420 times of maximum discharges, satisfy clinical requirements.



Numerous events review, trend and data storing.



With manual, automatic, power-on self-test function, ensure application at anytime.



80mm thermal printer, waveform is more clear and accurate.



Support 240 minutes of AED recording storage, the recording of each patient can be saved up to 60 minutes.

1 knob, mode can be switched among manual defibrillation, pacing and AED.

25 types, energy selections.

1s, electricity can be set up, rescue time is guaranteed fully.